
# See More. See It Now.

A deeper understanding through microscopy

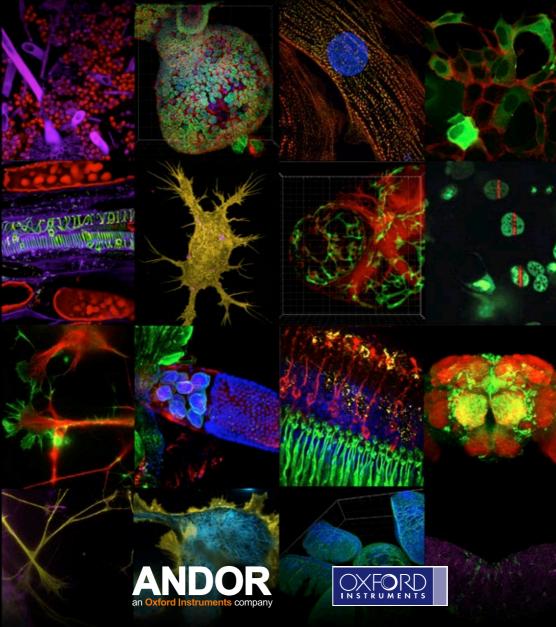

# **Dragonfly** - Confocal Imaging Platform

# It's more than confocal...

Andor's pioneering new Dragonfly confocal imaging platform offers multiple high-contrast fluorescence imaging techniques in a single device.

#### Instant confocal



Live confocal image (no line-by-line scanning)

- Up to 20x faster imaging for greater productivity
- Minimal photoxicity and photobleaching
- · Ideal for live and fixed sample imaging
- · High speeds for imaging fast dynamic events

#### Widefield imaging with accelerated deconvolution



Image samples not suited to other modalities

- The most effective technique for thin specimens
- Achieve sharp images with ClearView deconvolution
- · Fewer moving parts for faster wavelength switching

#### Simultaneous dual color TIRF



High contrast imaging at the membrane

- Match penetration depth of two labelled targets
  - Control penetration depth, critical angle, work in HiLo
  - Perfect for intra or extracellular imaging

#### Additional Features

Borealis Perfect Illumination Delivery™

Pinhole matching for low and high magnifications

Zoom illumination

3 selectable camera magnifications

Dedicated software control

Efficient user interface

Real-time 3D rendering

GPU accelerated deconvolution

#### **Key Benefits**

Extended imaging range into NIR

Broader choice of fluorophores

Deeper imaging

Best in class image quality

Outstanding field uniformity

Higher contrast images

Handles small and large samples

Ease of use and high productivity

Immediate experimental evaluation

#### **Key Applications**

Live and fixed cell imaging

Developmental biology

Cancer research

Organoids

Neuroscience

Membrane dynamics (e.g. exocytosis)

Single molecule localization imaging

#### Fusion Software - Acquire, Visualize, Analyze

Fusion is a brand new solution designed to meet the requirements of today's expectations for ease of use and immediate visual feedback with data review, whilst fulfilling tomorrow's aspirations for handling multi-modal imaging.

Fusion simplifies the control of the Dragonfly system, with its multiple imaging modes, to fluorophore and imaging mode selection in just a couple of mouse clicks. Once the sample is on the microscope you can control all hardware, including the x,y stage with its unique software joystick. Real-time 3D visualization provides a powerful insight to your experiment, and GPU accelerated deconvolution delivers enhanced clarity when required.

Saving files in Imaris format, Fusion permits easy transfer of the data into Imaris software for detailed multi-dimensional downstream analysis including solid surface rendering, measurement packages for cell biology, cell lineage, neuroscience and much more.

### Dragonfly - A multi-modal approach

- As an individual researcher you have the imaging tools to investigate from the whole organism down to a subcellular level, using multiple approaches to address a wide range of questions.
- If you run a core facility, you have a cost effective way to support multiple imaging techniques.

Find out more at andor.com/dragonfly





## Revolution DSD2

### The personal confocal imaging unit

Andor's Revolution DSD2 is an innovative hybrid of spinning disk technology and structured illumination. This unique approach is laser-free and delivers a budget friendly confocal solution to your laboratory, offering less dependency on laser based solutions that are often restricted to core facilities.

A simple device, which can even be added to an existing fluorescence microscope in your lab, the Revolution DSD2 will benefit your research by delivering confocal images as a routine technique in your work.

Whilst laser-free, the Revolution DSD2 can still achieve the optical sectioning and image quality you expect of a complex laser scanning confocal system, but with low maintenance costs. Furthermore, it does not need an expert to run it!

#### **Key Applications**

Immunofluorescence (3D structure and visualization)

Development e.g. C. elegans, Drosophila and Zebrafish


Embryology

Stem Cell and 3D cultures

Neuroscience

Cell and System Dynamics

Find out more at andor.com/dsd2

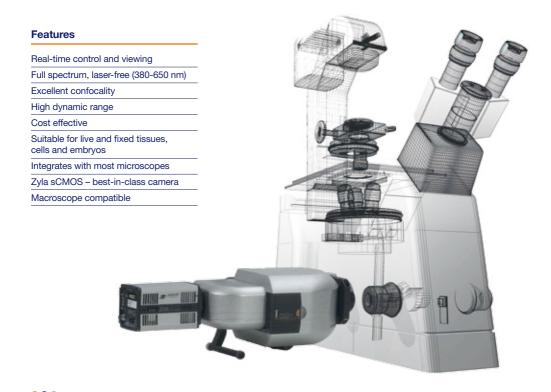

# Photostimulation Portfolio - Mosaic

### Simultaneous multi-region light targeting

The Mosaic active illumination system utilizes digital mirror device (DMD) technology to control the illumination field of a fluorescence microscope. Using a choice of illumination sources, Mosaic achieves real time and near diffraction limited resolution.

Unlike traditional galvoscanning systems, where pixels are addressed sequentially, Mosaic can simultaneously and precisely excite multiple regions of interest with complex geometries and allow simultaneous imaging.

#### **Key Applications**

Optophysiology/optogenetics

FRAP and PA

Studies of cell dynamics

Uncaging - signaling

Selective region imaging

Operating from 365-800 nm, Mosaic is unique yet flexible. Mosaic SDK offers access to software independent high speed pattern sequencing, ideal for applications such as optogenetics that mimic high speed cell signaling. Greyscaling is also available for detailed pattern illumination such as required in photolithography.

#### **Features**

Unlimited flexibility in shape and complexity of illumination mask

No scanning - simultaneous illumination of multiple regions of interest

Create complex pattern sequences with software independent recall

High speed pattern recall (up to 5,000 fps)

Multiple light sources available

Find out more at andor.com/ps





## Photostimulation Portfolio - MicroPoint

### Laser illumination and ablation

MicroPoint provides a flexible and field-proven tool for photostimulation. Supplied with a patented compact, pulsed nitrogen pumped tunable dye laser, it is capable of ablation, bleaching and uncaging over a wavelength range of 365 to 656 nm.

More than 20 wavelengths can be utilized with available dye resonator cells, and simultaneous imaging and photo-stimulation of the specimen can be achieved. MicroPoint is supplied with a UV-Vis imaging quality Epi illumination adapter for both current and previous generation microscopes from Leica, Nikon, Olympus and Zeiss.

There are two MicroPoint versions available:

- Manual simply centre in target and shoot a single diffraction limited spot. No software required
- Motorized software controlled x,y galvo mirrors steer the beam for region specified targeted light.

#### **Key Applications**

Ablation - Damage and repair
FRAP and PA - cell dynamics
Uncaging - signaling
DNA damage - cancer research
Semi-conductor marking
Circuit isolation




#### **Features**

Simultaneous laser delivery and image acquisition

Ablation, uncaging, activation and bleaching

365 - 656 nm - adapt to specific targets

Incremental control of pulse energy

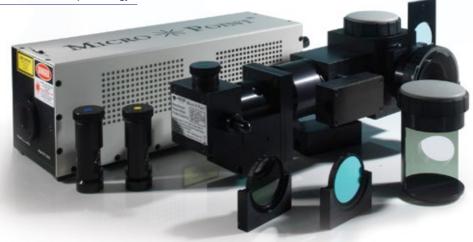

# Photostimulation Portfolio - FRAPPA

### Photo-bleach and activate

Andor's FRAPPA uses a dual galvanometer scan head to provide a computer-steered laser beam delivery system. By utilizing the ILE's range of lasers the FRAPPA provides unrivalled Fluorescence Recovery After Photo-bleach (FRAP) and Photo Activation (PA) flexibility.

In sharing the same lasers with imaging, you can bleach more efficiently close to the peak absorbance of the target fluorophore, and save money.

The FRAPPA provides diffraction limited performance and can be configured in-line with a spinning disk confocal or on other ports if so required. With high-speed galvo-driven mirrors, FRAPPA is capable of scan rates much faster than that of MicroPoint for photobleaching and photoactivation studies.

#### **Key Applications**

Photo-bleaching

Photoactivation

Vesicle tracking

Cell compartmentalization

DNA damage

Uncaging - signaling

Find out more at andor.com/ps

#### **Features**

All laser lines available for FRAP and PA actions

<10 ms switch over from Imaging to FRAP

Arbitrary multi-region scanning of points, rectangles and polygons

Integration with flexible protocols for 3D

FRAPPA localization and analysis



## Customer Support

Andor products are regularly used in critical applications and we can provide a variety of customer support services to maximize the return on your investment and ensure that your product continues to operate at its optimum performance.

Andor has customer support teams located across North America, Asia and Europe, allowing us to provide local technical assistance and advice. Requests for support can be made at any time by contacting our technical support team at andor.com/support.

Andor offers a variety of support under the following format:

- On-site product specialists can assist you with the installation and commissioning of your chosen product
- Training services can be provided on-site or remotely via the Internet
- A testing service to confirm the integrity and optimize the performance of existing equipment in the field is also available on request.

A range of extended warranty packages are available for Andor products giving you the flexibility to choose one appropriate for your needs. These warranties allow you to obtain additional levels of service and include both on-site and remote support options, and may be purchased on a multi-year basis allowing users to fix their support costs over the operating life cycle of the products.

#### **Head Office**

7 Millennium Way Springvale Business Park Belfast BT12 7AL Northern Ireland

Tel: +44 (0)28 9023 7126 Fax: +44 (0)28 9031 0792

#### **North America**

300 Baker Avenue Suite 150 Concord, MA 01742 USA

Tel: +1 860-290-9211 Fax: +1 860-290-9566

#### Japan

5F IS Building 3-32-42 Higashi-Shinagawa Tokyo 140-0002

Tel: +81-(0)3-6732-8968 Fax: +81-(0)3-6732-8939

#### China

Unit 1, Building A, Qing He Shun Shi Jia Ye Pioneer Park, No. 66 Zhu Fang Road, Haidian District, Beijing 100085 P. R. China

Tel: +86 (0)10-5129-4977 Fax: +86 (0)10-6445-5401

#### Find us on















